

CINCOM®
All For The Healthcare





 Cincom
 www.cincomhealth.com
 Email: service@cincomhealth.com

SHENZHEN CINCOM E-COMMERCE CO.,LTD
 Add: B806,11th Building ,Gongyuandadi Community,Longxiang Road,Longcheng Street,Longgang District,Shenzhen,China
 Web: www.cincomhealth.com
 E-mail: service@cincomhealth.com
 Made in China








CINCOM®
All For The Healthcare

Percussion Muscle Massager



Operating Instruction

Instruction

Thanks for purchasing our product. Please read the instructions carefully before use and keep it well for further reference. If you have any questions during using, please feel free to contact us through Amazon or mailing us by service@cincomhealth.com.

Safety Precautions

Warnings & Cautions

- Those with any of the following conditions or persons who are receiving medical treatment should consult the doctor before using the machine:
 - Using pacemaker or other medical devices which are susceptible to electrical interference;
 - Suffering from malignant tumors;
 - Suffering from cardiac diseases;
 - Having serious peripheral neuropathy dysfunction or sensory disturbance caused by diabetes;
 - Being unsuitable to do the massage because of traumas on the body;
- Keep it out of the reach of infants, children and people without the ability to use it independently.
- Do not immerse in water. Keep away from liquids or heat sources. Keep ventilation ports free from dust and debris.
- Do not remove screws or attempt to disassemble.
- Do not drop or misuse the Percussion Massage Gun.
- Do not use the machine continuously for more than one hour, stop using it until you have a rest for 30 minutes.
- Unplug the charger from the device after charging or before use.
- It can only be used on soft tissues of the body, which will not cause pain or discomfort. Do not use on any hard or bony area of the head or body.
- Fingers, hair, or any other part of the body should be kept away from the pole and back of the machine to avoid clipping.
 - Do not touch this part when the massage head is not installed and the device is on
 - Do not touch this part when the massage head is installed and the device is on


A面

10. Abrasions may occur during use no matter how much intensity you choose. Always check the area for treatment and stop at the first sign of pain or discomfort.

Feature

- Powerful Brushless High-Torque Motor with Quiet Glide Technology
- Rechargeable Lithium-Ion Battery (4-6 Hours of Use Per Charge)
- Multi-Speed Settings (4 Speed) for Optimal Pressure
- Cordless and Lightweight (1 Lbs)
- One-button-control
- Interchangeable Massage Heads (x4)

Product Diagram

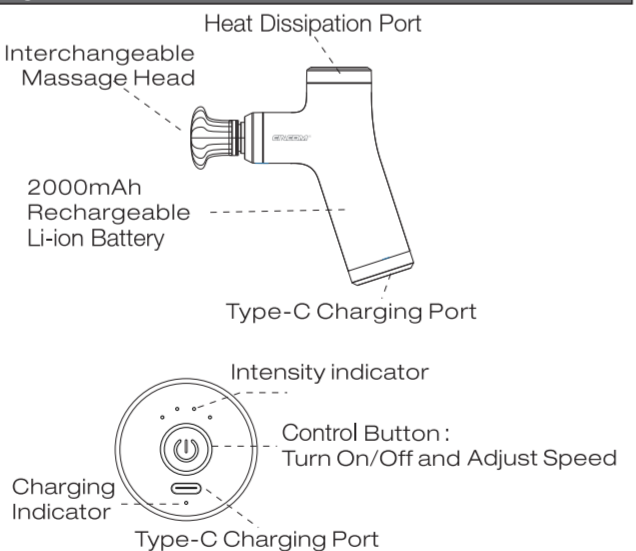

Operating

- Press the power button and hold for three seconds to turn on with level one. The blue lights indicate the speed level.
- Adjust vibration speed via power button. Maximum speed is level 4, minimum is level 1.
- Press the power button and hold for three seconds will turn off the device.
- Fully insert applicator by pressing in. Remove it by pulling out.

Notes: The device automatically turns off when using for 15 minutes. Press the power button again to continue using.

Charging Instruction

- For first charging, please make sure the machine is fully recharged, which is good for the battery.
- If the light of the power button is flashing, it means that the machine is out of power. Please stop using it and charge it.
- To recharge, connect the charging adapter to the device's charging port.

- When it is recharging, the indicator of power button is red and keeps flashing, and when it is fully recharged, the indicator light changes from red to blue. It can't be used during the charging.
- The average running time is 4+ hours depending on speed level and pressure applied during use.

Maintenance, Cleaning & Storage

- Use a slightly damp towel to clean device and a dry with a soft cloth.
- For storage or travel, Turn off the bottom power button and store in supplied portable box.

Package Included

- 1 x Massage Gun
- 1 x USB-C Charging Cable
- 1 x Operating Instruction
- 4 x Massage Heads

Specifications

Model	CM-052G
Name	Mini Percussion Muscle Massager
AC/DC Adapter	100-240V 50/60Hz, DC Output:5V 2A
Battery Capacity	2000Mah
Working Hour	4-6H
Charging Time	3-4H
No-load Speed	1500-2500 Rpm
Speed Levels	4

Contact Us

We provide 2-year warranty and lifelong after-sale service for all products in or store. If you have any problems while using, feel free to contact us via email.
Email: service@cincomhealth.com.

If you want to follow up our new products and latest promotional information, please visit our official website and Facebook account.

B面